

mediven® angio



Safety in treatment – down to the toes.

medi tip²:

Due to product features and the results of the clinical trial¹ we recommend mediven angio for:

- Women / men with chronic venous disease
- Concomitant mild to moderate peripheral arterial disease (PAD)
- Concomitant Diabetes mellitus







Indications³

- Chronic venous disease stages I bis III according to Widmer or C0s – C5 according to CEAP

Features


- Unique and soft merino-plush padding in the entire foot
- Natural merino wool for optimum thermoregulation and additional wearing comfort
- Extra long and flat seam clings to the toes and thus forms an anatomic tip shape for optimum pressure relief at the toes and their base joints
- Solid ribbed structure for an ideal balance of resting- and working pressure and as a visual guide for the correct application
- Micro-plush in the entire lower leg for a high wearing comfort
- 90°-heel ensures an even distribution of the knitted fabric in the heel and instep area

Washing instructions

Easy care    

Suitable for washing machine and dryer

medi tip:

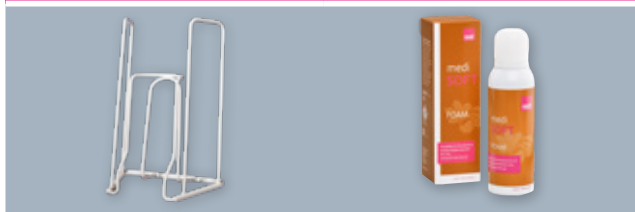

**medi Butler donning
made easier⁴**

**medi skin and textile care
offer support**

More information on the medi accessories from page 112

^{*} Product safety proven by clinical trial: The use of the medical compression stocking mediven angio in patients with chronic venous diseases and concomitant mild to moderate peripheral artery disease (PAD) and / or Diabetes mellitus is safe.¹

¹ Rother U et al. Safety of medical compression stockings in patients with diabetes mellitus or peripheral arterial disease. *BMJ Open Diab Res Care* 2020;8:e001316.

² Please note that this recommendation does not replace a doctor's diagnosis or professional advice from a specialist medical retailer.

³ Selecting the ideal compression stockings depends on the indication and the individual needs of the patient.

⁴ In general the following applies for donning and doffing aids: Rabe E et al. Medical compression therapy of the extremities with medical compression stockings (MCS), phlebological compression bandages (PCB), and medical adaptive compression systems (MAC): S2k guideline of the German Phlebology Society (DGP) in cooperation with the following professional associations: DDG, DGA, DGG, GDL, DGL, BVP. *Hautarzt* 2021;72(2):137-152.

mediven[®] angio


The solution for patients with chronic venous disease and concomitant PAD and / or diabetes mellitus.



Extra wide cuff

Secure hold and large area pressure distribution

- ✓ Prevention of constrictions

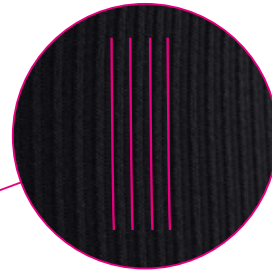

Solid ribbed structure

Easy-grip structure with high longitudinal elongation and fine plush material

- ✓ Easy to apply and comfortable to wear



Merino-Plush in the entire foot

Unique, soft padding from the tip of the toes up to above the ankle

- ✓ Reduction of pressure peaks in sensibility or skin disorders and optimum thermoregulation

Extra long, flat seam





Special seam which ends behind the base joints of the toes

- ✓ Optimum pressure relief of the toes and their base joints



mediven[®] angio

Safety in treatment – down to the toes.

mediven angio				
Styles	AD	AD	AD pte.	AD pte.
				
	calf small foot	calf large foot	calf petite small foot	calf petite large foot
CCL 1	3X40 Δ 11-7	3X40 Δ 21-7	3X47 Δ 11-7	3X47 Δ 21-7
CCL 2	3Y40 Δ 11-7	3Y40 Δ 21-7	3Y47 Δ 11-7	3Y47 Δ 21-7

Δ = Please add **1** for colour caramel and **5** for colour black
Please choose the size on the last digit of the article number from 1 to 7.

Standard styles	
Standard sizes	Size I – VII Length normal und petite
Made-to-measure	possible
Special models	not possible
Border	extra wide cuff
Foot part	closed toe small foot up to 25 cm $\hat{=}$ shoe sizes 36 – 40 large foot from 25 cm $\hat{=}$ shoe size 40 – 46

Material composition			
	Polyamide (PA)	Elastane (EL)	Merino wool
CCL 1	71 %	22 %	7 %
CCL 2	70 %	23 %	7 %



Material

Medium

Compression class

CCL 1, 2

Product benefits



Clima Comfort



Clima Fresh



Perfect Fit



Variety



Standard colours

S = Standard

M = Made-to-measure*

Black
S+M

Caramel
S+M

* Delivery period up to 10 working days.

Note

For more information please visit:
medi.biz/angio

Unique novelty

mediven® angio – the first and only medical compression stocking whose safety was scientifically evaluated and proven in patients with chronic venous disease and concomitant PAD and / or diabetes mellitus.*

Unique product features

mediven angio was particularly developed for patients with chronic venous disease and mild to moderate PAD and / or diabetes mellitus and optimally adapted to their needs.

Unique safety

In a clinical trial¹ with both types of patients, the safety of mediven angio was evaluated and scientifically proven:

- ✓ Stable blood circulation down to the smallest blood vessels
- ✓ No damage to the skin
- ✓ No constrictions
- ✓ No pressure points
- ✓ High wearing comfort

Confirmed for compression classes 1 and 2



* Product safety proven by clinical trial: The use of the medical compression stocking mediven angio in patients with chronic venous diseases and concomitant mild to moderate peripheral artery disease (PAD) and / or Diabetes mellitus is safe.¹

¹ Rother U et al. Safety of medical compression stockings in patients with diabetes mellitus or peripheral arterial disease. *BMJ Open Diab Res Care* 2020;8:e001316.